ITSU

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User Manual
IS0124 使用手册

User Manual

All information is correct at the time of print

Thank you for purchasing ITSU massager.

Please read the Safety Precautions carefully before using the foot massager. The Safety Precautions are divided into 'warnings' and 'attention' in terms of the degree of danger. This book should be kept available for future reference after reading.

CONTENTS

Safety Instructions	03
Warranty Claims	06
Product Parts	07
Product Features And Operation Illustrations	08
Product Operation Guide	09
Product Maintenance	11
Troubleshooting	12
Product Specifications	13
Contact Us	14

Safety Instructions

a) Important safety instructions

- To ensure a safe and proper operation, please read the instruction manual carefully before using this product.
- Use this product only as described in this manual.
- Always unplug this product from the electrical outlet immediately after use and before cleaning.
- Do not use this product if you have sensitive skin, neuropathy or poor circulation. Don't let children or people with mental disabilities use the heating function without proper supervision. Improper use may be dangerous.
- Do not operate this product under a blanket or when covered. Over-heating can occur and cause fire, electric shock or injury.
- Do not cover this product in order to avoid electrical faults or damages due to excessively high motor temperature.
- Stop the massage immediately if you feel any discomfort.
- Never operate this product if it has been dropped, has a damaged cord or plug; experiences a malfunction, or if it has been dropped in water. Please contact an authorized service centre to do the necessary check and repair.
- Do not use the components that are not recommended by ITSU.
- Do not insert any objects into any openings in the product.
- If you want to cut off the power supply, please make sure you have turned off all the operating controls, as well as remove the power plug from the socket.
- Please do not use this product as a medical device. If you need a professional diagnosis or treatment, please consult a professional medical person.

b) Operating Environment Precautions

- Do not use this massager in humid or extremely dusty environments to avoid malfunction.
- Do not use this massager at a temperature above 40 degree centigrade.
- Do not use this massager outdoors.
- Avoid placing or operating this massager in direct sunlight or near a heating appliance, stove.

Safety Instructions

c) User precautions

- Make sure that this product is far away from any hot surfaces.
- Please consult your doctor before use for those who have the following diseases:
 - Those who have heart disease, malignant tumor, back pain, spinal scoliosis, hypertension, osteoporosis, etc.
 - Those who are pregnant or having a menstrual period.
 - Those who have a physical discomfort, fever, a swelling or inflammation, bruises, allergies, etc.
 - Those who are warned by the doctor not to do any massage to avoid injury or physical discomforts such as blood clots (thrombosis), severe aneurysm, acute venous tumor, all kinds of skin inflammation, infectious skin disease (including subcutaneous tissue inflammation), etc.
 - Warnings: Do not use the massager together with the following medical instruments.
 - 1) Those who are using a pacemaker or other implantable medical equipment.
 - 2) Those who are using an artificial heart.
 - 3) Those who are wearing an electrocardiograph for similar medical equipment.

Otherwise it will cause malfunction for those medical equipment, putting the user's life in danger.

d) Precautions for power supply

- Do not turn on the switch when your hands are wet.
- Do not immerse in or spray water or liquid on the control power switch.

e) Cord/Plug Precautions

- Do not connect or remove the plug with wet hands. Otherwise you will suffer an electric shock or a product malfunction.
- Do not pull the wires when you pull out the power plug. Please hold the plug itself.
- As the default operating time is 20 minutes, all the functions will stop working after 15 minutes. Please unplug this product from the electrical outlet immediately after use.
- Keep the cord away from hot surfaces.
- Do not place the power plug under the product or any heavy items.

Safety Instructions

- Do not wind the power plug around the product lest it should be damaged and cause an electric shock or fire.
- Do not operate this product if it has a damaged cord or plug.
- Do not use this product if the power socket is loose.
- Do not carry this product by its supply cord or use the cord as a handle.
- When the power plug is damaged, please call the service centre or an authorized technical person for the necessary check and maintenance.



Attention

- Do not fall asleep while massager is on.
- Do not use this product while you are drunk or physically uncomfortable.
- Don't stand on the product or put a heavy weight on it.
- Do not use this product for more than 30 minutes.
- Do not use this product if the cloth cover is worn out.
- Do not use this product within an hour before or after a meal.
- When you are using the product at a low room temperature, do not suddenly use the heat function. You can adjust the room temperature to normal room temperature, before starting the heat function.
- When moving this product from a low temperature storage place to a warmer setting, it may cause malfunction due to water vapor gathered in the components. Set the product in a normal room temperature for one hour before using it.
- Do not use this product continuously if you find any components damaged or the product in an abnormal condition.
- When current interruption occurs during your massage, try to avoid further damage due to sudden power restoration.
- Stop using the product if there is thunder and lightning.



Warranty Claims

1. Warranty is not transferable.
2. This warranty is effective only if the product is purchased from ITSU.
3. Proof of purchase (original receipt) is required for all warranty repairs.
4. All implied warranties, including but not limited to these implied warranties of fitness and merchantability, are limited to the duration stated in the warranty policy.
5. Damage due to misuse, improper treatment and unauthorized modification and repair are not covered by this warranty.
6. Warranty does not cover accessories and attachments which do not belong to the massage unit.
7. Warranty is not effective to rental, business, commercial, institutional or other non-residential users.
8. All services covered by this warranty must be approved by ITSU and repaired by authorized technicians only.
9. If replacement parts for defective materials are not available, ITSU reserves the right to make substitution in lieu of repair or replacement.

Product Parts

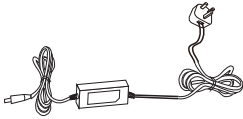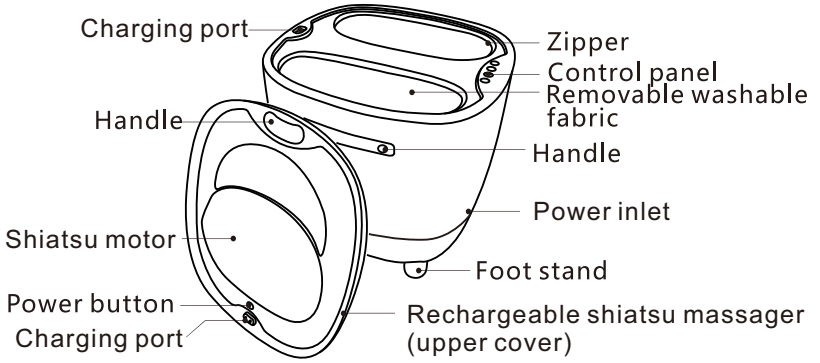The power adapter

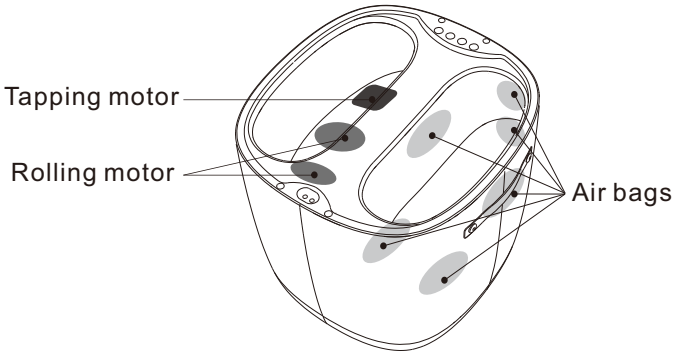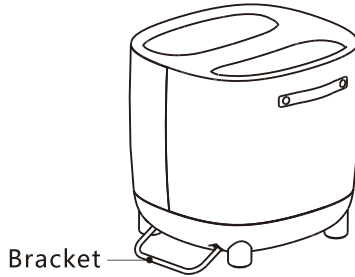Product Features And Operaton Illustrations


Preparation before use

- ▶ Please put the product on the flat floor, and then remove top cover.

Usage method

Step 1: wire connection

- 1 Connect power adaptor to the power line plug-hole located on the product as shown on the right.
- 2 Plug the power plug into the AC 100-240~V power outlet.

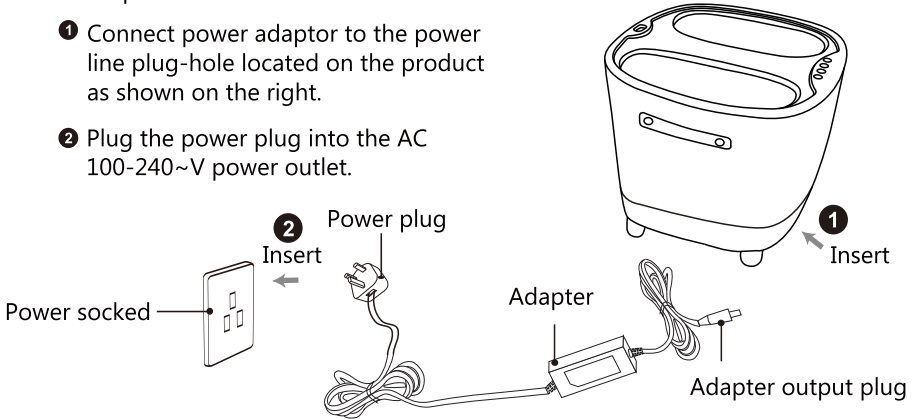Caution

- ▶ Before plug, please check and make sure the power cord is in good condition.
- ▶ Do not use the wet hand to pull out the power cord in order to avoid any electric shock.

Step 2:

- ▶ Please choose a chair with desired height and then put your feet into the product. Turn on the power button on the control panel to enter the massage state. Take off the top cover cushion and place it on a chair or sofa cushion. Press and hold the power button to start kneading massage. You can also hold your stomach or other part of the massage. (See right picture 1) According to the use method 2, the upper cover cushion can be closed to charge.

Picture 1

Picture 2

Step 3:

- ▶ The upper cover can be used as a stool, (see picture 2).

After use

- ▶ Turn off the power button on the control panel.
- ▶ Remove the plug from the power supply socket when not in use.

Caution

- ▶ Please do not touch power connect points after cover heated in case of a hot.

Product Operation Guide

ON/OFF

Press to turn ON.
Press again to turn OFF.
Default setting timer is 15 minutes.

Auto mode button

Press to enable AUTO program 1 to 5.
Program 1: Airbags massage @ blue light
Program 2: Rolling & tapping massage @ red light
Program 3: Series of soft massage @ green light
Program 4: Series of moderate massage @ purple light
Program 5: Series of powerful massage @ white light

Heat button

Press to turn ON.
Press again to turn OFF.

Intensity button

Press to switch message INTENSITY
Weak(green)>Middle(blue)>Strong(red)

Operation Instruction for Shiatsu Cover

Charge for 2.5 hours, can be used for 2 hours continuously.

Battery charging indicator:

Charging: red light illuminate

Fully charged: green light illuminated

Operation method:

Press on hold the POWER button for 2 seconds to start massage with heat

Press POWER button again to change massage direction

Press POWER button again to turn off heat function

To turn off the product, press POWER button for 2 seconds again.

*Massage direction changes automatically every one minute

*The massager turns off automatically after 15 minutes.

Caution

- ▶ Please do not use the product by standing posture or put any heavy object on the product to avoid any damage of the product.
- ▶ Please do not stand or seat on the product.
- ▶ When moving the product ,please do not pull the power cord or drag the product directly to avoid any damage.
- ▶ Do not try the product on a no-load basis to avoid wear and tear.

- Each massage program is 15 minutes.Product will stop once massage is completed.

Product Operation Guide

a) Instructions before use


Please store this product in a dry place, do not expose the product to sunlight
put the carpet under the product to minimize the working noise.

b) How to use it

This product uses only the AC 220V-240V, 50/60Hz power supply.

Make sure your indoor voltage fits this specification
parameter.

- Connect the AC power cord to the socket.

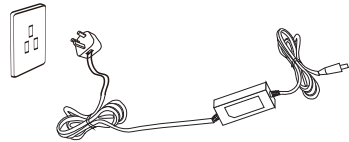Do not pull the plug
with wet hands.

Product Maintenance

a) Storage

- Please store it in the packing carton, and put in a safe, dry, well-ventilated, and shaded place.
- Do not twist the power plug.
- Do not hang the power wire and pull the product.
- Avoid sharp items from destroying the surface of the product.
- Please don't store the machine under direct sunlight or in places of high temperature.

b) Clean

- Before cleaning, make sure that the power cord is unplugged from the electrical outlet and that the machine is completely cool.
- Use a damp cloth to clean the product. Ensure the product is completely dry before reuse.
- The detachable fabric cloth can be washed in a normal washing way.

*Attention

- Make sure not to use benzene or thinner and other corrosive detergents to clean the product.

c) Maintenance

- Keep the product clean and store in a dry environment.
- When the product is not in use for a long time, please put a dustproof cover on the product.

Troubleshooting

If you meet the following problems, please referring to the following solutions. If the product is still malfunctioning or doesn't work, please send the product to our service centre or an authorized technician for the necessary checks and repair.

Conditions	Possible reason	Solutions
It doesn't work	<ol style="list-style-type: none">1. Plug isn't connected correctly to the socket.2. The power switch at the back is off.3. Switch key on control panel isn't on.	<ol style="list-style-type: none">1. Push the plug firmly into the socket.2. Switch on the power switch at the back of the product.3. Switch on the on/off key on the panel.
Stop working suddenly	<ol style="list-style-type: none">1. The power switch is off.2. Massage time is over 15 minutes.	<ol style="list-style-type: none">1. Check that the power plug is connected well.2. Restart the power supply switch.

Product Specifications

Power supply	AC 220V-240V, 50/60Hz
Power consumption	60W
Timer	15 minutes
Dimension(LxWxH)	430x400x480mm
Weight	12Kg
Dimension of packing box (LxWxH)	465x450x515mm
Weight (including the packing box)	13.2Kg

Contact Us

ITSU World Sdn Bhd
13 Lorong Teknologi B. Nouvelle Industry Park.
Taman Sains SG1 Kota Damansara. PJU 5.
47810 Petaling Jaya
Malaysia
Tel: +603 78903366

ITSU World(HK) Limited
Workshop D1,28/F. TML Tower
No. 3 Hoi Shing Road, Tsuen Wan, N.T.
Tel: +852 3611 2228

PT Healthy World
Lodan Center Block O No 20
Jalan Lodan Raya No 2
Jakarta Utara 14430
Tel: +62 21-6983 1123

User Manual

All information is correct at the time of print

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