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User Manual  
IS2002 使用説明書

User Manual  
All information is correct at the time of print

# CONTENTS

Safety Instructions	03
Warranty Claims	05
Main Unit & Accessories	06
Product Operation Guide	07
Product Maintenance	09
Troubleshooting	10
Contact us	12

# Safety Instructions

## **i. Usage environment**

- Do not use the appliance in a wet condition or dusty environment as this may result in malfunction or electric shock.
- Place the appliance in a safe location.
- Do not expose the appliance to a heater or direct sunlight.
- This product is only intended for personal use only. Commercial use is prohibited.

## **ii. Health precautions**

If you are experiencing health complications or conditions as stated below, please consult your doctor before use:

- Suffering from heart disease.
- Pregnancy.
- Any health concerns.

## **iii. When the product should not be used**

- If the use exceeds 10 minutes at one time.
- This unit should not be used by children or invalids without adult supervision.
- Before sleeping.
- If the fabric cover is torn.
- If water is accidentally spilled onto the appliance.
- If the user feels abnormal or severe pain during the massage.
- If the user detects any malfunction or other abnormal conditions during operation.
- If there is lightning.

## **iv. Caution and warning**

- Unplug the appliance from the outlet when not in use.
- Do not use any attachments that are not recommended by ITSU.
- Do not operate this appliance if it has a damaged cord or plug. Please return the unit to ITSU Service Center for examination and repairs.
- Do not insert your hands or any foreign matters into the massager's vents and gaps.
- Handle heated surfaces carefully. The unattended use of the heated surfaces by children or invalids may be dangerous.
- This product is only intended for personal use only. Commercial use is prohibited.

## Safety Instructions

- Do not cover the appliance during operation. As this may result in malfunction or fire.
- Do not crush the appliance and avoid any sharp folds.
- The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.

### V. Note

- This appliance can be used by children at the age of 8 or above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge under guidance or if they have been informed of the risks involved and instructions concerning use of the appliance in a safe way.
- Children shall not play on the appliance.
- Cleaning and maintenance shall not be performed by children without supervision.
- To ensure safe and proper use of the appliance, do not operate it before reading this user manual carefully first.
- It is normal to hear the noises during use, because they are resulted from vibration of product.
- For people under medication or with a medical condition, please consult your doctor before use. Please do not use the product if you are not feeling well.
- If the appliance is used in a cold room, do not increase the room temperature abruptly. It is advised that the temperature is gradually increased to a normal level.
- If the appliance is stored in a cold place and being taken into a warm environment, please wait one hour before use to get the optimum massaging level.
- Do not remove the cover of the appliance. Mishandling the internal components may result in malfunction or electric shock.
- Do not connect or disconnect the power plug from the electrical outlet with wet hands.
- Do not use the appliance if the electrical outlet is loose.
- The appliance is only to be used with power supply provided.
- Children shall not play with the appliance.
- Cleaning and user maintenance shall not be made by children without supervision.


## Warranty Claims

1. Warranty is not transferable.
2. This warranty is effective only if the product is purchased from ITSU .
3. Proof of purchase (original receipt) is required for all warranty repairs.
4. All implied warranties, including but not limited to these implied warranties of fitness and merchantability, are limited to the duration stated in the warranty policy.
5. Damage due to misuse, improper treatment and unauthorized modification and repair are not covered by this warranty.
6. Warranty does not cover accessories and attachments which do not belong to the massage unit.
7. Warranty is not effective to rental, business, commercial, institutional or other non-residential users.
8. All services covered by this warranty must be approved by ITSU and repaired by authorized technicians only.
9. If replacement parts for defective materials are not available, ITSU reserves the right to make substitution in lieu of repair or replacement.

# Main Unit & Accessories

- A) Massager
- B) Power adapter
- C) Car adaptor
- D) Wireless remote

## Controller Diagram

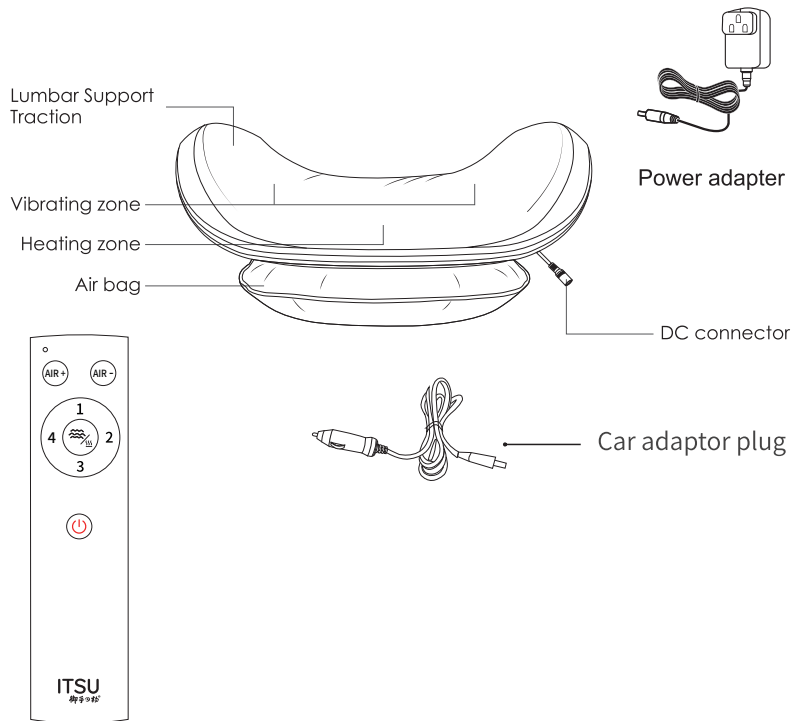
# Product Operation Guide

- i .People with back diseases should consult a doctor before using the product. When using the product for the first time, we recommend you manually inflate or deflate the air bag, and 50% of full inflation is enough.
- ii .Heating function: It takes about 1 minute for the product to warm up.
- iii. When using the product for the first time, please do not inflate the air bag excessively to avoid hurting the waist.

1. Lie down slowly. After lying and adjusting the position, begin the massage and control the working modes via the controller.

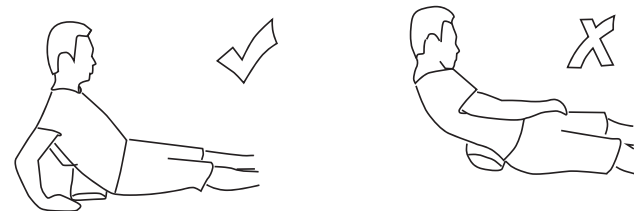2. The product will shut off automatically after 10 minutes of massage, or you can stop the massage in advance by pressing the power button. After the massage has finished and the air bag has been completely deflated, continue to lie for 1-3 minutes before getting up slowly with support of arms.

## ii. Operation instruction

### a) Preparation before use

- Check if there is any damage to the appliance.
- Place the appliance on the part which needs to be massaged, before ensuring that the appliance is properly placed.

### b) Usage method

1. Connect the DC head of the power adapter to the DC connector of the appliance, and plug the other side into the power base.

● Default timer is 10 minutes

● Key operation: Each effective button press will come with 'beep' sound.

 Press to ON/OFF the product. Notes: Vibration, heat function can be enabled or disabled in auto program

 Press to enable Auto 1, a gentle massage program.

 Press to enable Auto 2, a moderate massage program.

 Press to enable Auto 3, a stronger massage program + vibration.

 Press to enable Auto 4, a stronger massage program + vibration + heat therapy.

 Manual mode: Press to enable vibration 1, vibration 2, vibration 3, OFF.

Auto mode: Press to add or disable vibration in auto program.

Once Power ON, long press for 2 seconds to enable or disable heat.

 Once Power ON, long press for 2 seconds to control air inflation time.

 Once Power ON, long press for 2 seconds to control air deflation time.

# Product Maintenance

## TO STORE

- Place the massager in this box or in a safe, dry and cool place. Avoid contact with sharp edges or pointed objects that might cut or puncture the surface.
- Please wind up and store the power cord under the appliance bottom frame. Do not hang the unit by the power cord.

## TO CLEAN

- Avoid immersing the unit in any liquid during cleaning.
- Unplug the unit and allow it to cool before cleaning.
- Clean only with a soft, slightly damp sponge.
- Do not allow water or any other liquid to go inside the unit.
- Avoid using abrasive cleaners, brushes, gasoline, kerosene, glass/ furniture polish or paint thinner to clean the product.
- Do not attempt to repair the product. For servicing, please contact ITSU.

## CAUTION

- Avoid tearing on the appliance surface.
- Only gentle force should be exerted against the unit in order to eliminate any risk of injury.
- Beware when put on light color PU/fabric seat, some material may have color dyeing
- Beware when put printing document on product, some material may have color dyeing



# Troubleshooting

This massager is controlled by an integrated chip and is susceptible to power surges from +10% to -15%. Due to this, the massager might display the following:

Conditions	Possible reason	Solutions
NO Inoperation	Poor power outlet	1. Check the power outlet for electrical output 2. Check whether the outlet plug is in good contact.
Power socket is good but not working	fault	Please unplug the plug and contact the seller for repair.

# Product Specifications

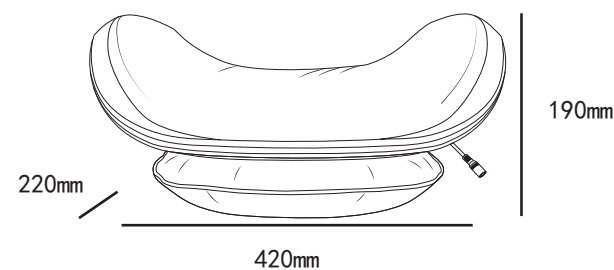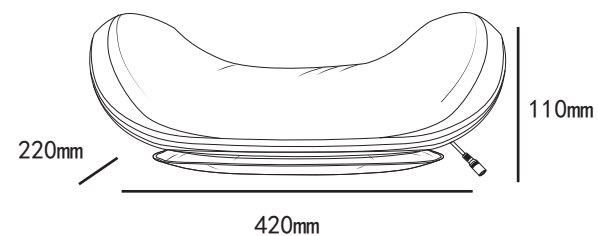Model number of the product	IS2002
Power source	100-240V AC
Power frequency	50/60Hz
Power consumption	12W
Voltage	12VDC, 1A
Timer	15 minutes
Dimension (Length/width/height)	420*220*110 (mm)

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ITSU World(HK) Limited