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Rev.1.0

User Manual
IS0602 使用説明書

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User Manual

All information is correct at the time of printing

Thank you for purchasing ITSU product.

Please read the Safety Precautions carefully before using the foot massager. The Safety Precautions are divided into 'warnings' and 'attention' in terms of the degree of danger. This book should be kept available for future reference after reading.

Contents

Safety Instructions	2
Warranty Claim	4
Product Parts	5
Product Operation Guide	7
How to use	10
Product Maintenance	12
Common Problems & Solutions	16
Product Specifications	18
Contact Us	19

Safety Instructions

Before operating the unit, please read all instructions carefully.

Make sure that the product is connected to a dedicated electric outlet for 220-240Vac, 50 Hz/60Hz.

Electricity Safety

- Do not handle the power plug or electric outlet with wet hands.
- Do not use a damaged power cord or plug.
- Plug in the power plug properly to operate the product.
- Do not attempt to disassemble or repair any electronic parts of the product.
- Do not spill water or liquids on the product.
- Do not use the product in a damp or humid environment.

Operation Safety

- Do not use the product near heating areas or on top of a heating surface.
- Do not use the product if you are experiencing health complications or conditions as stated below:
 - Suffering from heart diseases, malignant tumors, back pain, spine disorders, osteoporosis, high blood pressure and knee joint degeneration.
 - Pregnant or menstruating
 - Feeling unwell, fever, having inflammation, bruises and sensory impairment.

Others

- When there is a power outage, carefully unplug the product.
- If the product functions abnormally, immediately turn off the power, disconnect the power plug and contact our service centre.
- If you feel unwell or experience any discomfort while using the product, stop using it immediately.
- Close supervision is necessary when the product is near children, infants, or disabled persons.
- Do not insert your hand or foot into any openings of the product.
- The product is intended for use of one person each time.

Safety Instructions

Warning

- Do not fall asleep while using the product.
- Do not use the unit after drinking alcohol.
- Suggested usage time is 20 minutes, please rest for 30min before use it again.
- if you feel unwell or sick, Please stop immediately.
- Please follow the instructions accordingly while using the unit.
- While the product is in use, do not insert your hands into the space on the two sides of the product.
- Children should be supervised to ensure that they do not play with the product.

Use caution when handling the product, as some parts may contain sharp edges or protrusions that could cause injury if not observed.

Correct Disposal of this product

This marking indicates that this product should not be disposed with other household wastes. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources.


Caution: read this user manual carefully before use.



Warranty Claim

1. Warranty is not transferable.
2. This warranty is effective only if the product is purchased from ITSU or any authorized distributors.
3. Proof of purchase (original receipt) is required for all warranty repairs.
4. All implied warranties, including but not limited to these implied warranties of fitness and merchantability, are limited to the duration stated in the warranty policy.
5. Damage due to misuse, improper treatment and unauthorized modification and repairs are not covered by this warranty.
6. Warranty does not cover accessories and add-ons, which do not belong to the product.
7. Warranty is not effective to rental, business, commercial, institutional, or other non-residential users.
8. All services covered by this warranty must be approved by ITSU and repaired by authorized technicians only.
9. If replacement parts for defective materials are not available, ITSU reserves the right to make substitutions in lieu of repair or replacement.

Product Parts

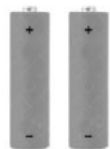Product Parts

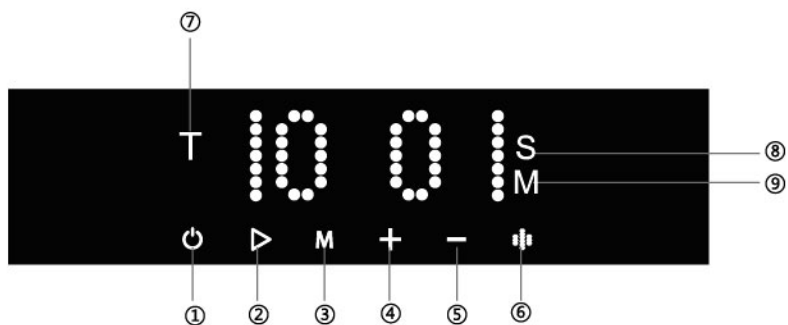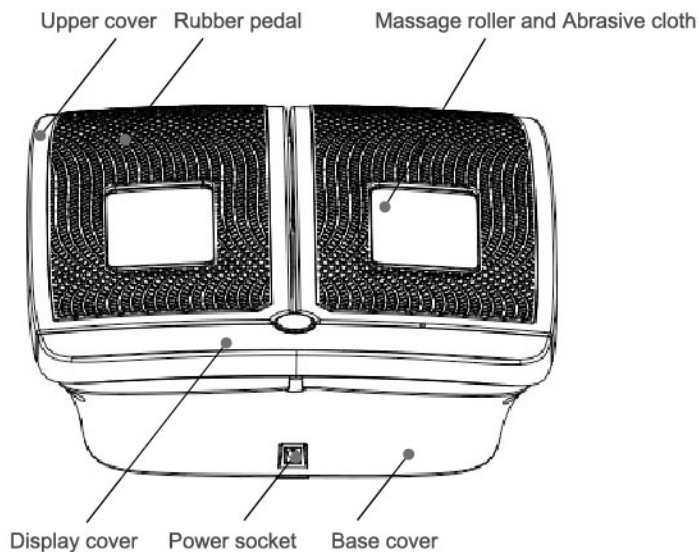
- | | | |
|------------------|----------------|---|
| ① Power Button | ② Start/stop | ③ Movement (with roller or without roller) |
| ④ Increase speed | ⑤ Reduce speed | ⑥ Auto mode |
| ⑦ Time | ⑧ Speed | ⑨ Display M |

Product Operation Guide

Please check before use

The product is for 220-240Vac, 50 Hz/60Hz.

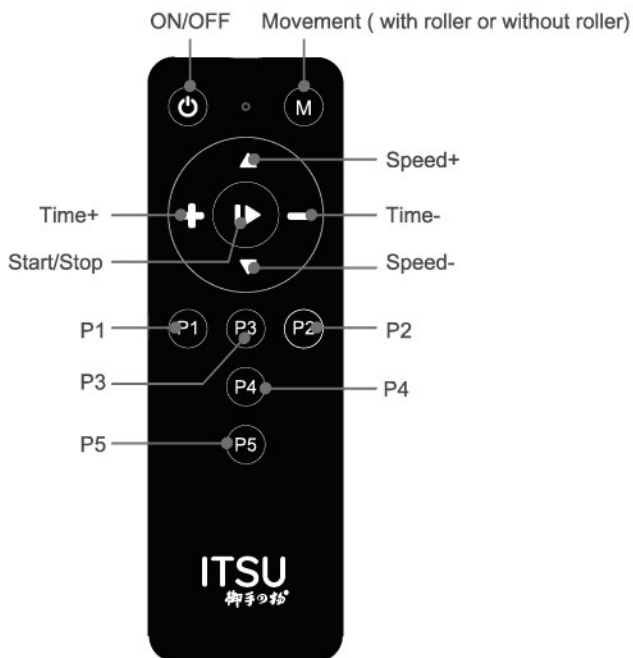1. Power :keep power supply energized.
2. start/stop : control start and stop of the machine running state.
3. Movement : control the two major movement of the machine, press the button to choose with roller or without roller (just up and down).
4. Speed increase: Manual 1-20. Press this button. To increase the speed when the machine is operating.
5. Speed decrease: Manual 1-20. Press this button.To decrease the speed when the machine is operating.
6. Time increase:Maximum 20minutes, press this button to increase time. (the device's default setting is 10 minutes).When the machine is operating.
7. Time decrease: Maximum 20minutes, press this button to decrease time.(the device's default setting is 10 minutes).When the machine is operating.
8. P1,P2,P3,P4,P5 : Auto program Select Auto program by default and press the start/stop button to start.

Product Operation Guide


Speed range 1-20 levels

A=movement without Roller


B=movement with roller

preset program 10 mins each, after press started the time will countdown.

P1

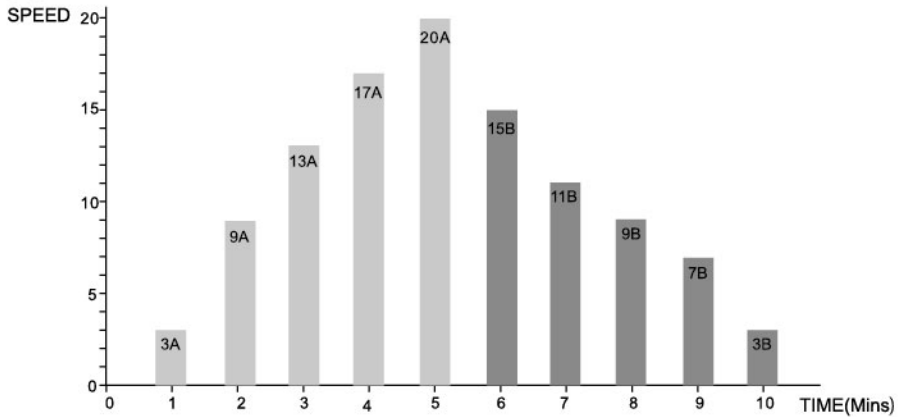
P2

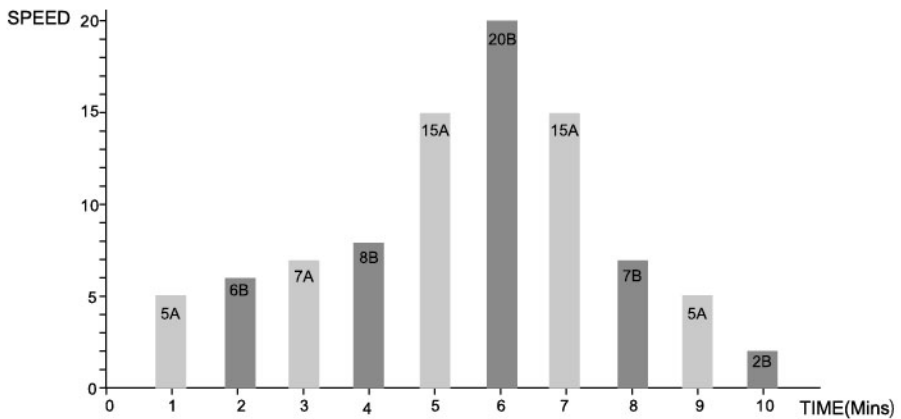
Product Operation Guide


P3

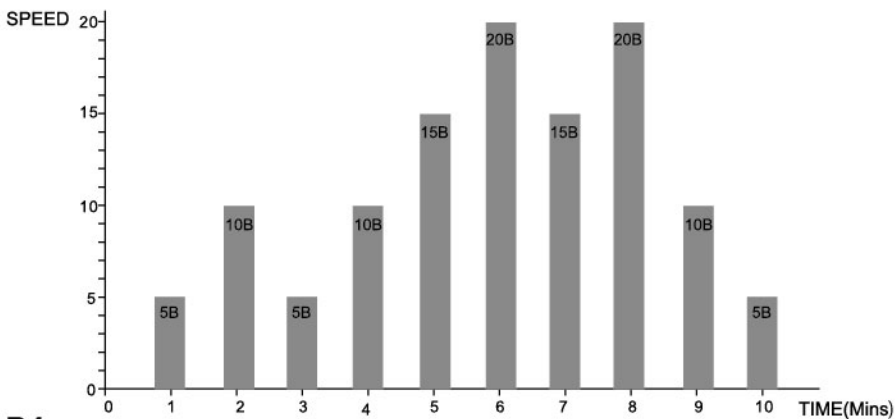P4

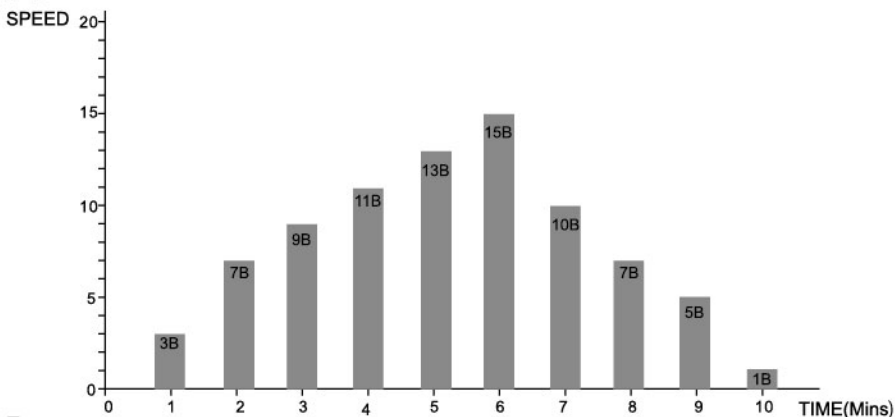P5

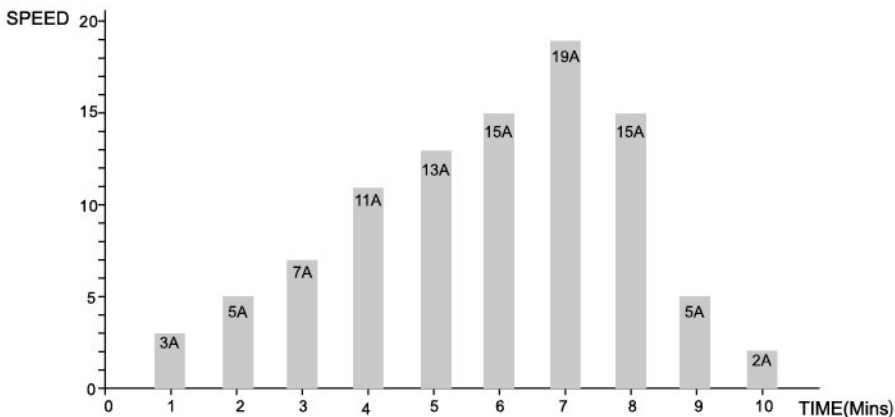How to use

place your feet and finally, set your speed. There you go, now you can get all the benefits of this passive gym wherever you go.

·Its movements are precise, continuous and with vibrations, this will keep your legs moving continuously, improving your blood circulation, helping you prevent leg cramps and back pain while also decreasing vascular problems and leg inflammation.

·Improve your health: Creates constant leg movement to mimic the natural muscle activity of walking in the convenience of your own home or office.

·Improve your figure: Our leg exercise machine also helps you improve your figure by fighting with cellulitis and helping you to prevent those awful varicose veins. It also strengthens your muscles everytime you use it. Just sit back, relax and with just 20 minutes a day you will see results.

1. Foot vibration massage

Under the office desk

Sit on the couch

Lateral varus of foot

Tiptoe cock

Use on tiptoe


How to use

2. LEG vibration massage


Calf/belly use

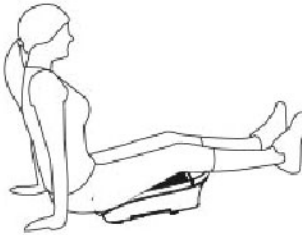Use of inner thigh


Back and Waist massage

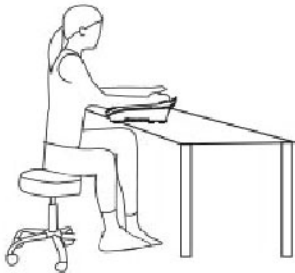Arm vibration massage

Hand vibration massage

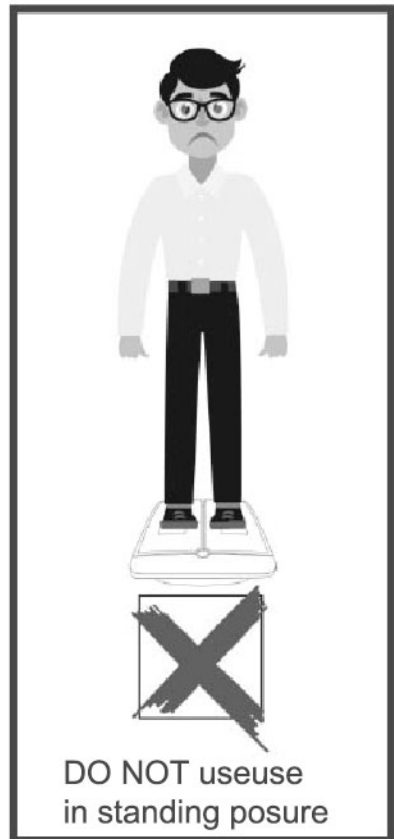Product Maintenance

- If the product was stored in a cold place for a period of time, it is suggested to wait one hour before usage.
- If the product is going to be placed in a special environment, such as salty or sulfurous area, please contact our service center for protection to prevent malfunction.
- Always unplug the unit before cleaning it.
- Please do not use thinner, benzene or alcohol to clean the unit.
- Use dry cloth to thoroughly clean the product, do not wipe with wet cloth.
- Do not place the product near heaters, or exposed to heat resources or direct sunlight.
- If you may not use the product for a period of time, it is suggested to unplug the unit and store it well under a clean and dry environment.
- Do not eat or drink on the product, to prevent malfunction.
- Storage environment :Temperature-10°C~55°C, Humidity RH0~85%

AUTO	1	2	3	4	5	6	7	8	9	10
P1	3A	9A	13A	17A	20A	15B	11B	9B	7B	3B
P2	5A	6B	7A	8B	15A	20B	15A	7B	5A	2B
P3	5B	10B	5B	10B	15B	20B	15B	20B	10B	5B
P4	3B	7B	9B	11B	13B	15B	10B	7B	5B	1B
P5	3A	5A	7A	11A	13A	15A	19A	15A	5A	2A

Speed range 1-20 levels

A=movement without Roller

B=movement with roller

preset program 10 mins each, after press started the time will countdown.

Product Maintenance

※When cleaning, be sure to turn off the power and unplug the power cord from the socket.

Wipe dirt with a dry, soft cloth.

※ if the dirt is serious, please put the diluted neutral detergent in the rough cloth, wring it out and then wipe it. After that, please wipe it with a dry cloth.

Be sure to clean the remote with a dry cloth, as cleaning the remote with a wet cloth may cause malfunction.

Do not clean the equipment with water. Do not use diluent, benzene, alcohol, remover and other organic solvents to clean the host. Otherwise it may fade, corrode and crack.

Product Maintenance

Please store this product in the following places

		Temperature 10 ~35 °C Humidity: 20 ~80%
Horizontal and stable place	Ventilated area	Above temperature and humidity location

Please do not store this product in the following places

The device is made of sophisticated mechanical and electronic components. Avoid installing the machine in the following places, which may cause a failure.

		
Places exposed to direct sunlight	Places with a lot of dust	Places where temperature changes rapidly
		
Places where the humidity changes drastically	Places with fire	Places that are easily wet by water
		
Places with volatile substances	Places close to cooling and heating equipment	Places with vibration

Caution

Do not use it on infants or people who cannot express their intentions, and do not keep it within easy reach. Otherwise it could lead to an accident or poor health.

Product Maintenance

- If the product was stored in a cold place for a period of time, it is suggested to wait one hour before usage.
- If the product is going to be placed in a special environment, such as salty or sulfurous area, please contact our service center for protection to prevent malfunction.

Always unplug the unit before cleaning it.

Please do not use thinner, benzene or alcohol to clean the unit.

Use dry cloth to thoroughly clean the product, do not wipe with wet cloth.

Do not place the product near heaters, or exposed to heat resources or direct sunlight.

- If you may not use the product for a period of time, it is suggested to unplug the unit and store it well under a clean and dry environment.
- Do not eat or drink on the product, to prevent malfunction.
- Storage environment :Temperature-10°C~55°C, Humidity RH0~85%

Common Problems & Solutions:


Q1	Does not work when power is turned on	Q5	Can it be used in the bathroom?
A	Check that the power cord is connected to both the host and socket. Check whether the main power switch is on. Check the operation panel or remote control is not working. * if only the remote control does not work, please replace the battery, because the remote control battery may be dead.	A	This product is not waterproof. Please do not use it in the bathroom as it may cause trouble.
		Q6	Can it be used outdoors?
		A	Don't use it. Please do not use it outdoors as it may cause trouble.
Q2	Suddenly stopped	Q7	Can I use it during pregnancy or right after giving birth?
A	Check that the power cord is connected to both the host and socket.	A	Please do not use. For postpartum and lactating persons, consult a physician who frequently visits the clinic.
Q3	Make noise	Q8	Is there a weight limit?
A	Be sure to install it in a flat and stable place. The structure of the product produces some operational noise, so it is not an anomaly in this case. If the noise is heard inside the product, foreign body may be mixed in due to damage or other reasons. Please contact the customer service window on page 00.	A	The load capacity of the product is 100kg. Do not use if the load capacity is exceeded.
			Remote control does not respond
		A	Please replace the battery of the remote control. Other reasons may be the remote location of the product and remote control, not suitable for display panel. Please aim the remote control at the display panel and operate within a distance of 1.5m.
Q4	Can I use it every day?	Q10	I lost my remote
A	No problem but only use it three times a day (30 minutes) at most. Overuse can lead to accidents, such as injuries.	A	Consult your dealer.
Q11	I don't know where to attach the band	Q16	Stuck, unable to operate
A	On the back of the main frame there is a connecting plate for the upper body band. Please refer to page 24 for details.	A	Replace the battery of the remote control and check the operation of this product again. In addition, the product is designed to stop automatically 10 minutes after use, either automatically or manually.
Q12	I want to attach the ground terminal, but there is no mounting hole at home	Q17	Strong current, feel pain
A	After consulting with a dealer or electrical construction shop, we recommend ground installation work (class D < type 3 > installation work).	A	The output setting may be too high. Try adjusting the output using the remote control. In addition, EMS mode is used only once a day.
Q13	Weak stimulation in EMS mode		
A	Output Settings may be low. Adjust the output with the remote control.		
Q14	I feel the stimulus is different on the left and right		
A	If only one foot is wet, or if you use a moisturizer, the current experience will be different. Please do not use this product until you have thoroughly wiped your feet. Also, even if only one leg has a scratch or cut, the body will feel different. Do not use this product until the wound is completely healed.		
Q15	Warning sound does not stop * Manual (EMS mode)		
A	Due to dry skin, the current may not flow easily. Please use again when the contact surface is wet		

Common Problems & Solutions:

S/N	Problems	Causes	Solutions
1	The Product can not be started	The plug is not in place	Insert the plug
		The power switch is off	Turn switch on, the indicator
2	Sound like phew or Da-da is heard during working	Normal sound generated by the Motor or from the abrasion Between mechanical units	No additional treatment is required
3	Operate the controller, but the Product no action	The electricity of batteries is worn out	Replace batteries
		The controller is too far away from the Product	The effective distance of infrared signal is 2.5M
4	Sudden stop	Default working time is up	Restart it after 10 minutes'interval
5	Can not operation	Too many commands sent to the Product and they can not be responded in time	Power off and then restart
6	Display "ER01" or "ER02"	The motor's working current is over	Disconnect the power supply, then Insert the plug and restart or contact with the after-seller

Product Specifications

Voltage	220 – 240Vac, 50/60Hz
Rated power	50W
Model	IS0602
Size	W380mm × D320mm × H135mm
Weight	Approx 4.5kg
Remote control power supply	2 pcs AAA batteries
Material	Body: ABS resin, rubber Remote control: ABS resin Special mat: TPE

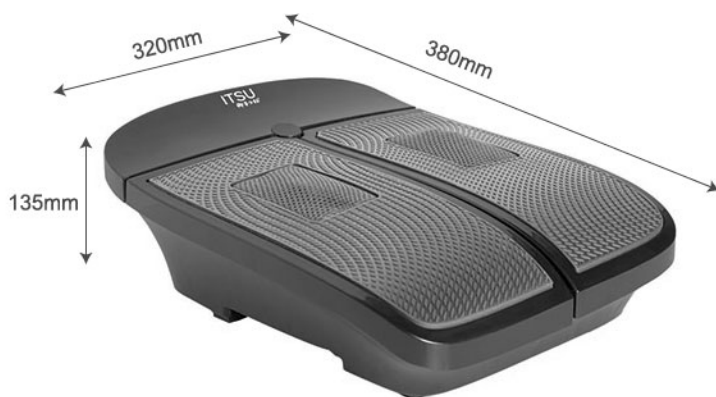Contact Us

ITSU WORLD SDN. BHD. (reg. no. 1137036-K)

Unit BO03-B-13A-1,
Menara 3A, KL Eco City,
No.3, Jalan Bangsar,
59200 Kuala Lumpur.

Tel: 1700-81-4878

ITSU World(HK) Limited

Workshop D1,28/F. TML Tower
No. 3 Hoi Shing Road, Tsuen Wan, N.T.

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